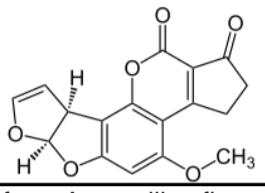




Specification Sheet

Antigen Aflatoxin B1-BSA		
Product Name	AFT B1-BSA	
Catalog Number	B94-Ag1	
Lot Number		
Description		Aflatoxin B1 Bovine Serum Albumin Conjugate (AFT B1-BSA conjugate) from <i>Aspergillus flavus</i>
Product Profile	Link position specification	AFT B1-BSA Conjugate
	Appearance	Lyophilized powder, faint yellow
	solubility	Soluble in water up to 2 mg/ml, (in most application, 10 mM PBS pH 7.4, with 0.1% NaN_3)
	Hapten/protein	~10:1
	Purification	chromatography purified
	Concentration	
	Purity	>95%
	Other	
	Application	Suitable for use in ELISA or Lateral flow immunoassay. Other applications have not been tested but use in such assays should not be excluded.
Performance review	<ul style="list-style-type: none"> Recommended pair for competitive inhibition immunoassay: Coating: B94-Ag1; Labeling: A94-Ab1, A94-Ab2 	
Store at	Refrigerate at 2-8°C	
Precautions and Disclaimer	This product contains NaN_3 , which has been classified as Harmful, in European Directive 67/548/EEC in the concentration range of 0.1–1.0%. When disposing of this reagent through lead or copper plumbing, flush with copious volumes of water to prevent NaN_3 build-up in drains.	
Note	Centrifuge before opening to ensure complete recovery of vial contents	
FOR RESEARCH OR FURTHER MANUFACTURING USE ONLY, NOT FOR USE IN DIAGNOSTIC PROCEDURES		